
beper

| | | |
|----|--|---------|
| IT | MASSAGGIATORE A PERCUSSIONE RICARICABILE | pag. 2 |
| EN | RECHARGEABLE PERCUSSION MASSAGER | pag. 6 |
| FR | MASSEUR À PERCUSSION RECHARGEABLE | pag. 10 |
| DE | WIEDERAUFLADBARES PERKUSSIONSMASSEGERÄT | pag. 15 |
| ES | MASAJEADOR DE PERCUSIÓN RECARGABLE | pag. 19 |
| GR | ΕΠΑΝΑΦΟΡΤΙΖΟΜΕΝΗ ΣΥΣΚΕΥΗ ΜΑΣΑΖ ΜΕ ΚΡΟΥΣΗ | pag. 24 |
| RO | APARAT DE MASAJ PRIN PERCUȚIE REÎNCĂRCABIL | pag. 28 |
| CZ | DOBÍJECÍ MASÁŽNÍ PŘÍSTROJ | pag. 32 |
| NL | OPLAADBAAR MASSAGEAPPARAAT | pag. 36 |
| LV | UZLĀDĒJAMA PERKUSIJAS MASĀŽAS IERĪCE | pag. 40 |

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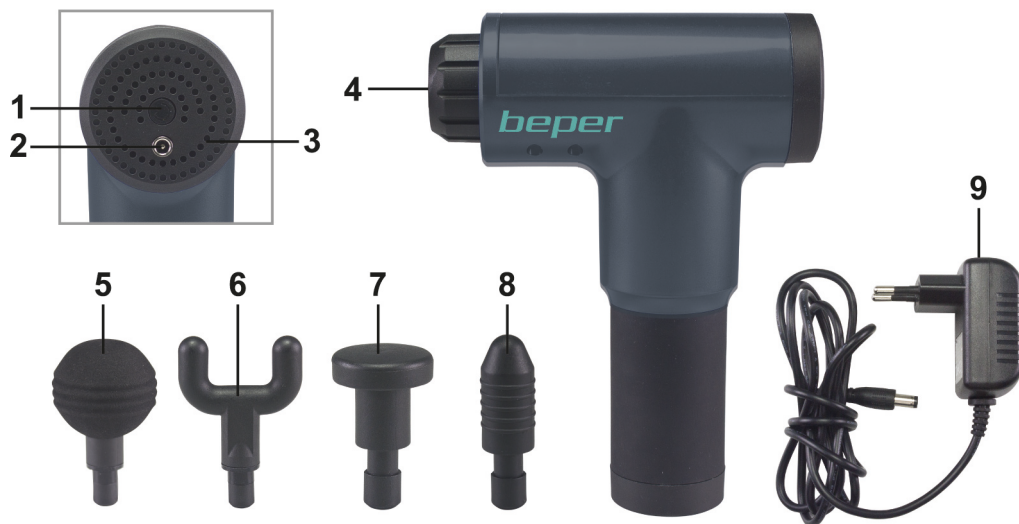assistenza@beper.com

General safety advice

Read the operating instructions carefully before first use of the appliance.

Read the following instructions carefully, since they give you useful safety information about installation, use and maintenance and help you avoid mishaps and possible accidents.

Remove the packing and make sure that the appliance is intact, with special attention to the supply cable. The elements of the packing (plastic bags, polystyrene, etc) do not have to be accessible to the children, in order to preserve them from danger sources.

It is recommended not to throw these packing material in the domestic waste, but to deliver them to the appropriate station of collection or destruction refusals, asking eventually information to the manager of the sanitation department of your own municipality.

This appliance is to be destined to its special purpose only.

Any other use is to be considered as improper and thus dangerous. The manufacturer cannot be responsible for any damages deriving from improper, wrong or incautious use.

To protect against electric shock, do not immerse any parts of the appliance into water or other fluids.

Do not expose the appliance to atmospheric agents such as rain, wind, snow, hail.

Do not use the appliance near explosives, highly flammable materials, gases, burning flames, heaters.

Do not allow children, disabled people or every people without experience or technical ability to use the appliance if they are not properly watched over. Close supervision is necessary when any appliance used by or near children.

Make sure that children do not play with this appliance.

If you decide not to use the appliance any longer, make it unserviceable.

Also make any dangerous parts unserviceable, especially for children who could play with the appliance.

Warnings for use

Do not immerse device in water. Keep device away from liquids and/or heat sources.

Keep ventilation holes clean and free from dust, dirt, etc. to avoid overheating of motor.

Do not remove components, screws, etc., and do not attempt to disassemble/repair the device.

Do not use the device for more than 1 hour continuously. It is advisable to turn the device off after 30 minutes of continuous use, to allow the motor to cool down, and then turn the device back on.

Make sure to always disconnect the charging cable before use.

This device cannot be used by children, or persons with reduced physical, sensitive or intellectual capabilities, or with lack of experience and knowledge, without supervision or if they are not given instructions and warnings on the correct and safe use of the device and the potential risks related to its use.

Do not use the device on open wounds or post-operative bruises; before use, consult your physician in case of any disease or physical condition that requires special attention (e.g. wounds, muscular trauma, pregnancy, etc.).

This device is not a medical device, and it is not designed to diagnose, treat or prevent any disease or medical condition.

While using, make sure to keep beard, hair, etc. away from moving parts of the device, to avoid them from getting trapped.

Do not use this device to treat areas above the Adam's apple, in the proximity of the cervical vertebrae (C3 and C4), on the head, or near the genitalia.

Do not use accessories other than those provided with the device and included in the package.

Only use the adapter provided with the device; use of different adapters can lead to damage of the device, and cause potential hazards to user's health and safety.

Product description Fig. A

1. ON/OFF and intensity adjustment button
2. Connector for charging cable
3. Ventilation holes
4. Connection hole for massaging accessories
5. Ball accessory
6. U-shape accessory
7. Flat accessory
8. Bullet accessory
9. Adapter and charging cable

Charging of device

Before first use, the device must be fully charged.

Connect the adapter to the device, by inserting the cable inside the connection hole found on the device, and then connecting it to the wall socket.

On the device, right below the connection hole, a LED light turns on (red color), while the light indicator on the adapter turns red; when charging is completed, the light on the adapter turns green: battery is now fully charged. Full charging requires approx. 2,5 hours, and allows up to 1,5 hours of use (battery life may vary based on intensity level selected and pressure applied during use).

Installing of massaging accessories

Choose the desired massaging accessory based on the area to be treated, and according to your personal needs.

The device comes with 4 different interchangeable accessories:

Ball accessory: for extended areas of the body, such as pecs, quadriceps, lumbar area, neck and foot plant

U-shape accessory: for a targeted and deep massage of the areas adjacent to the spine and for the scapular area

Bullet accessory: ideal for small areas, such as hands and feet, for the treatment of tense muscle groups and for the stimulation of trigger points

Flat accessory: for a wide and extensive massage, suitable for all areas of the body

Once the accessory has been selected, insert its connector inside the connection hole found on the front of the device, by applying a gentle pressure. Make sure that the connector of the accessory is completely inserted inside the connection hole. Gently pull outwards to remove the accessory.

⚠ Important: make sure the device is turned off before inserting/removing a massaging accessory.

Use

Always turn the device on before placing it on the area of the body to be treated.

Press the ON/OFF button once; the device turns on, working at the lowest massaging intensity level.

By repeatedly pressing the ON/OFF button, it is possible to increase intensity, from level 1 (the lowest) to level 6 (the highest). A red LED indicator will turn on above the ON/OFF button, indicating the intensity level selected. To turn the device off, remove it from the treated area, and press the ON/OFF button 7 times (or until the device stops working).

It is advisable to start the treatment with the lowest intensity level, and then to gradually increase level, to avoid any possible bruise or trauma.

Gently move the device around, applying a gentle pressure on the area to be treated, making sure not to stop on the same spot for more than 60 seconds.

Use the device only on dry and clean skin, free from clothes or other objects that might hinder the movement of the massaging accessory.

Do not use this device to treat areas above the Adam's apple, in the proximity of the cervical vertebrae (C3 and C4), on the head, or near the genitalia.

Immediately stop using the device in case of intense pain or discomfort.

It is advisable to turn the device off after 30 minutes of continuous use, to allow the motor to cool down, and then turn the device back on.

Cleaning and maintenance

Make sure that the device is turned off, and adapter is unplugged from device and wall socket, before proceeding with any cleaning and/or maintenance operations.

Clean the outer surface of the device with a soft, slightly damp cloth; wipe thoroughly before storing.

⚠ Warning: do not immerse device in water. Keep device away from liquids and/or heat sources.

Technical data

2000 mAh lithium battery

Adapter:

Input 100-240V 50/60Hz 0.5A

Output 8.4V ==== 1.0A

For any improvement reasons, Beper reserves the right to modify or improve the product without any notice.

The European directive 2011/65/EU on Waste Electrical and Electronic Equipment (WEEE), requires that old household electrical appliances must not be disposed of in the normal unsorted municipal waste stream. Old appliances must be collected separately in order to optimize the recovery and recycling of the materials they contain, and reduce the impact on human health and the environment.

The crossed out "wheeled bin" symbol on the product reminds you of your obligation, that when you dispose of the appliance, it must be separately collected.

GUARANTEE CERTIFICATE

This appliance has been checked in the factory. From the date of original purchase a 24 month guarantee applies to material and production defects. The purchase receipt and the guarantee certificate must be submitted together in case of claim to guarantee.

The guarantee is valid only with a guaranty certificate and proof of purchase (fiscal receipt) indicating the date of purchase and the model of the appliance

For any technical assistance, please contact directly the seller or our head office in order to preserve the appliance's efficiency and NOT TO void the guarantee. Any intervention on this appliance by non-authorized persons will automatically void the guarantee.

GUARANTEE CONDITIONS

If the appliance shows defects as a result of faulty material and/or production during the period of guarantee, we guarantee repair free of charge on condition that:

- The appliance has been used properly and for the purpose for which it has been intended.
- The appliance has not been tampered with, otherwise it is not maintainable.
- The purchase receipt shall be presented.
- The appliance showing a fair wear and tear shall not be covered by this guarantee.

Therefore, any part that could be accidentally broken or having visible signs of use in consumable products (such as lamps, batteries, heating elements...) the aesthetic parts are excluded from the guarantee, and any defect whatsoever resulting from non-respect of the rules for use, negligence in usage and/or maintenance of the appliance, carelessness, wrong or improper installation, damage during transport and any other damage not attributable the supplier.

For every defect that could not be repaired within the guarantee period, the appliance will be replaced free of charge.

In any case, if the part to be replaced for defect, breakage or malfunction is an accessory and/or a detachable part of the product, Beper reserves the right to replace only the very part in question and not the entire product

Contact your distributor in your country or after sales department Beper.

E-mail assistenza@beper.com which will forward your enquiries to your distributor.

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